

Latest
+
Greatest



InSite® In-Facility Medication Packaging and Dispensing System

- + Minimizes waste**
- + Transforms workflow efficiency**
- + Reduces pharmacy labor**
- + Shortens time to administration**

«The InSite System definitely saves us time on our med-pass. It saves us about 50% than previously so we have more time to spend with residents and give them quality care.»¹

Transforming Medication Management for Long-Term Care

The InSite System is changing the way medications are dispensed in long-term care facilities. This in-facility medication packaging and delivery technology enables secure, automated dispensing on location at long-term care facilities. InSite improves patient safety, reduces medication waste and saves nursing time by delivering pharmacist-approved, on-demand medication 24/7/365.

With InSite, the medication dispensing process is automated and accurate. Nurses can quickly package pharmacist-approved medications before each med pass, or on-demand to facilitate new residents, First Doses, STAT orders and PRNs.

SOURCE

1 — Long-Term Care Nurse, Milner Community Healthcare

Advantages

Minimizes Waste – With on-demand packaging, changes in prescribed medications or patient residency will no longer produce unused blister packs that need to be destroyed or returned. This reduction in waste processing delivers an average of \$4,000 savings per month per facility.

Transforms Workflow Efficiency – When contrasting in-facility medication packaging to the traditional blister card model, InSite significantly reduces the number of steps and minimizes the opportunity for human error.

Increases Revenue – Short cycle dispensing fees enable long-term care pharmacies to receive an average of \$1,500 per facility each month in additional dispensing fees for each facility that adopts in-facility dispensing.

Reduces Pharmacy Labor – On average, a long-term care pharmacy will save \$12,000 per year in pharmacy technician hours filling canisters instead of blister cards.

Shortens Time To Administration – InSite dramatically reduces the time needed to prepare for scheduled medication passes, freeing nurses to spend more time with their residents.

Benefits

- Eliminates blister pack-related storage and inventory management concerns.
- Dramatically reduces the time needed to prepare for scheduled medication passes.
- Eliminates the need to perform time-consuming narcotics counts.
- Works with existing information systems, basing its automated inventory replenishment requests on actual medication usage.
- Reduces the number of medication deliveries and expensive special courier runs required.



InSite How it works

With InSite, medication orders are transmitted to the pharmacist for clinical review. After approval, the medications are immediately available at the long-term care facility and dispensed from the InSite System with multi-dose, patient-specific packaging.

- Before delivery, each canister is filled, security-sealed, verified and signed by the pharmacist.
- Medications are dispensed only as needed and for a specific patient.
- Changes in prescribed medications or patient residency will no longer produce unused blister packs that need to be destroyed or returned.
- InSite provides inventory and dispensing information for down-to-the-dose accountability at each facility.



1

Order submitted to the pharmacy via phone, electronically or fax



Pharmacist performs initial prescription verification (PV1)

2


Order is sent electronically to the InSite System



InSite System

3

Nurse utilizes InSite System to package/dispense unit or multi-dose patient-specific medication

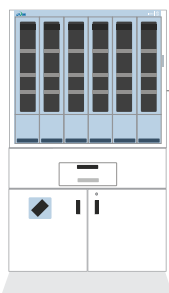

4

Nurse easily removes medication from pouch




5

Medication is administered to the resident



6

InSite System keeps record of dispense

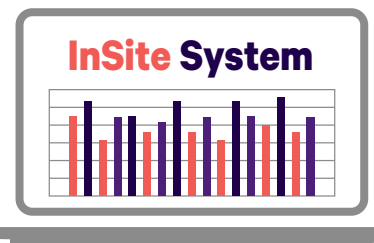

InSite® In-Facility Medication Packaging and Dispensing System

The InSite System stores up to 300 different medication types in a single medication dispensing unit, providing automated, on-demand compliance packaging for oral solid medications prior to administration.

1 Medication Security

Restricts access to stored medication with password protection and individual locking storage drawers. Staff access to medication is documented and logged.

2 High-Speed Packaging

Packages up to 60 doses per minute. Rapid dispensing for each cycle minimizes the wait time for staff.

3 Fully Front Loaded Unit

Canisters are easily stocked directly from the front of the packager. The required amount of workspace is minimized since there is no need for routine access to the sides of the packager.

4 Efficiency

The average InSite canister holds about 300 doses and can be filled and processed in roughly the same time it takes to create a single 30-day blister pack.

5 Swappable Lower Packaging Unit

Redundancy for critical components is assured with the ability to swap the packaging component on-site to minimize downtime.

6 Safety

Supports an Automatic Canister Recognition System (ACRS) for added convenience and packaging accuracy.



6 Inventory Management

Allows automatic inventory monitoring as well as lot and expiration tracking.

7 Medication Return Bin

Nursing staff can eliminate medication waste by securely storing returns in the electronically controlled and monitored return bin.

8 User Friendly

Color-coded guides simplify the loading of ribbon cartridges to eliminate misalignments and printing issues.

9 Smart Design

Overall footprint will fit through most standard doors without modification for easy installation. The floor brace ensures front stability while the front drawers are open.



Partnering with Swisslog Healthcare

We're transforming medication management for long-term care facilities and pharmacies.

At Swisslog Healthcare, we strive to lead change for better care. At the core of this vision is a focus on improving workflows and reducing the time clinicians spend doing routine and repetitive tasks—enabling more time to care for patients and residents. In addition to long-term care automation, our solutions

and services extend across the continuum of care, including transport, medication and supply chain automation for hospitals and health systems.

Headquartered in Buchs, Switzerland, Swisslog Healthcare has been delivering best-in-class automation solutions for more than 100 years. We are a global organization with 2,500 team members in more than 20 countries and clients in more than 50 countries.



CONTACT

Swisslog Healthcare
healthcare.us@swisslog-healthcare.com
800.764.0300
swisslog-healthcare.com/insite

InSite is a registered trademark of Swisslog AG. Swisslog Healthcare systems may be covered by one or more patents. See swisslog-healthcare.com/patents for details.
© 2022 Swisslog Healthcare EN NAM 102022 ISS-102